ApexRp 0.1 CFM Power Distribution Box

The R0xP Power distribution box powers the 0.1 CFM ApexRp without requiring an extra cable, providing a convenient solution. It serves as a compact junction box enabling power and data transfer

in a single package.



RJ45 Connector to RPC

Pin 4 (Data +) from the Particle counter will connect to terminal 1 of the junction box

• Pin 5 (Data -) from the Particle counter will connect to terminal 2 of the junction box

 Power wire from the power cable will connect to terminal 3 of the junction box

Ground wire from the power cable will connect to terminal 4 of the junction box

12345678

RJ45 Pinout

T-568B

1. White Orange

3. White Green

2. Orange


5. White Blue

7. White Brown

6. Green

RJ45 Connector to Control Cabinet/ Gateway .



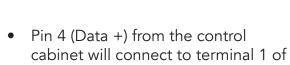

RPC

Power and Data Cable ~1 foot

in Length

• Pin 5 (Data -) from the control cabinet will connect to terminal 2 of the iunction box

the junction box

- Pin 7 from the control cabinet will connect to terminal 3 of the junction
- Pin 8 from the control cabinet will connect to terminal 4 of the junction box



Lighthouse Worldwide Solutions reserves the right to change specifications without notice.

