ActiveCountR


VIABLE IMPACTOR SYSTEM



IGHTHOUSE designed the ActiveCountR system to meet the industry need for continuous trouble-free operation and flexible running parameters. The ActiveCountR utilizes an external vacuum source which can be an independant pump system, house vacuum or a shared vacuum source with non-viable particle counters.


The sampling head is desinged to be remote mouted within the critical sampling environment. The sampling head comes in a bottom mount and side mount design. The remote sampling switch allows the operator to start, stop or pause sampling without having to intervene into the critical area being sampled. Sampling can also be controlled via software.


FEATURES

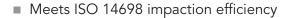

- 1 CFM flow rate
- Local or remote start/stop control via software
- External vacuum source for clean operation
- Stainless steel enclosure
- LED status indicators for power, flow and sampling
- Auto power fail recovery
- 316L stainless steel sampling head
- Designed for reliability
- 2 year warranty



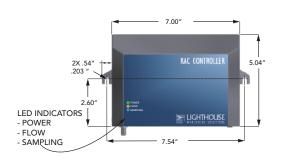

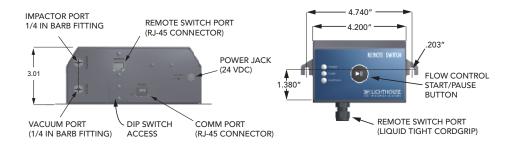

TECHNICAL DATASHEET

ActiveCountR Viable Impactor System

Features		
Flow rate	1 CFM.	
Sampling Methodology	Impaction Principle.	
Communication Protocol	RS-485 Modbus.	
LED Indicators	Power, Flow, Sampling.	
User Interface	Manual Start/Stop; Control Software (LMS Pharma, LMS Pro)	
Communication Connectors	Switch: Liquid-tite cord grip; Controller: RJ45	
Enclosure	Stainless Steel	
Sample Outlet Connection	1/4 " Barb	
Vacuum Requirements	16 in HG	
Power	24 VDC	
Dimensions Impactor	Diameter: 4.41"(11.20 cm), height 2.1" (5.33 cm)	
Dimensions Switch	$5.8"$ (w) x 2.7 (h) x $1.0"$ (d) $[15 \times 6.9 \times 2.5 \text{ cm}]$	
Dimensions Controller	W 8.0" (w) \times 5.0 (h) (20.32 cm), H 5.0" (12.7 cm), D 3.0" (7.62 cm)	
Weight Impactor	Impactor: 4.4 lbs (2.00 kg)	
Weight Switch	Switch 0.59 lbs (0.27 kg)	
Weight Controller	Controller: 2.6 lbs (1.18 kg)	

Environmental Conditions

 $\begin{array}{ll} \textbf{Operating} & 50^{\circ} \text{F to } 104^{\circ} \, \text{F} \, (10^{\circ} \, \text{C to } 40^{\circ} \, \text{C}) / 20\% \, \text{to } 95\% \, \text{non-condensing} \\ \textbf{Storage} & 14^{\circ} \, \text{F to } 122^{\circ} \, \text{F} \, (-10^{\circ} \, \text{C to } 50^{\circ} \, \text{C}) \, / \, \text{Up to } 98\% \, \text{non-condensing} \\ \end{array}$

Accessories

Included Installation guide on CD; power supply, sample tubing

Optional Printed installation guide; LMS Express, Express RT software, LMS Pharm software, LMS Pro software

Distributed By:



Lighthouse Worldwide Solutions reserves the right to change specifications without notice.

